User Manual

6-in-1 Massager

User Manual

6-in-1 Massager


Please read this user manual carefully to fully understand the usage of the 6 in 1 Massager. Keep the manual for future reference since it contains important information.

WARNING

1. Must pull out the power cord and connecting line when cleaning, maintenance and remote product.
2. Must stop working When exception occurs (burning taste).
3. Do not use wet hand to pull the power cord.
4. Do not use in humid place as the bathroom.
5. Do not spray coating, volatile oil and pesticides.
6. Do not use broken adaptor, if soft wire power or plug damaged, must use special soft line to replace.
7. Must clear the adaptor dust regularly.
8. Do not put the product in the sun exposure long-time.
9. Do not modify the product casually.
10. Do not let children near when working, dangerous.
11. Should cut off adaptor when going out or sleeping.
12. There should be people when working long-time.
13. Do not input combustible and do not use in oil.
14. Do not wash the product with water.

STRUCTURE

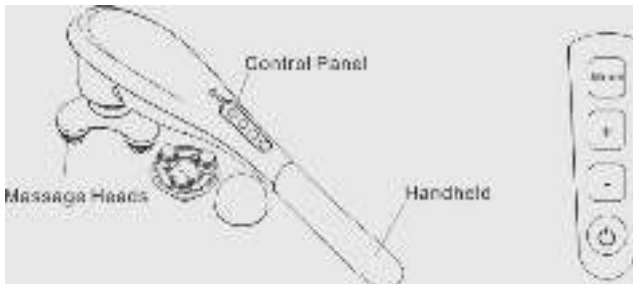CHARACTERISTICS

- Cabinet model and moderate massage strength
- Changeable massage heads, can massage body's special part
- Can be used in all parts of the body, eliminate pain
- Can be used for abdominal fat push, body beauty
- Solve constipation problem

Tripod massage head

Simulated hammering and finger pressing, Wide massage area, Suitable for massaging the whole body

U-shaped massage head

Simulate multi-point massage, Even strength and soft, Suitable for massaging neck, leg and other curved parts

Fan shaped massage head

Simulate multi-point massage, For Shoulders, Leg and other body parts

Round massage head

Traditional anthropomorphic massage, Suitable for massage all parts of the body

Needle massage head

Simulated acupuncture and massage, Suitable for stimulation and tapping Hand, sole and calf

Press the massage head

Simulate the feeling of pressing, Suitable for massage Specific parts

INTRODUCTION

- Insert the plug into the socket, press the power switch, 4 modes and 3 speed can be adjustable.
- Massaging 1~2 times every day, 5~15mins every time, can be health care.
- Finishing working or not in use temporarily, please off the switch, then pull out the power plug.

CLEANING & MAINTANCE

1. Do not pull the cord with hand, and do not let the product impending and shaking off, or the power cord core will disconnect, and not work normally.
2. Prevent water or other corrosive liquid enters into, the fault or damage caused.
3. Please use neutral cleaner and soft cloth to wipe, corrosive cleaner prohibited

TECHNICAL SPECIFICATIONS:

Rated Voltage: DC5V 1A
Rated Wattage: 20W
Vibration Partial Rate: 2.5mm
Vibration Frequency: 3200rpm~0

Imported/Distributed by The Investment Company Europe B.V. exclusively developed for sales in the Netherlands, Belgium, and Luxembourg.

PATENT PENDING

The artwork and design of this package are protected by European copyright law and may not be reproduced, distributed, displayed, published or used for any purpose without prior written permission. It is not permitted to alter, or reproduce any trademarks or copyright notice from this package.

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